



A device for disconnecting IV, dialysis and feeding connections

Provides users extra grip and leverage, so connections come apart with less force.

- Straight bar design gives added strength.
- Soft, non-slip jaws provide a secure grip on the connection without damaging the line.
- Three jaw sizes fit a wide range of IV, dialysis and feeding tube connections.

Indications for use

To limit force on hand and wrist, use Yewtwist on all IV disconnections

- For routine use on every disconnection, and on overtightened, difficult to disconnect or adhered fittings.
- Intended to assist with the disconnection of luer lock fittings, including fittings used for medical purposes, by providing additional grip and leverage.
- Do not use hemostats, clamps or other tools to disconnect.


Approved use

Single patient use only

For the highest level of aseptic practice, dedicate one Yewtwist to each patient for the full course of their therapy.

- Keep it at the bedside or with the patient's own supplies. Do not keep it in a shared drawer, cart or pocket.
- Write the patient's name on the label of the resealable bag.
- Clean before and after each use with hospital-grade wipes.
- Discard when the patient's therapy ends, or replace sooner if worn.
- Do not share between patients.

Directions for use

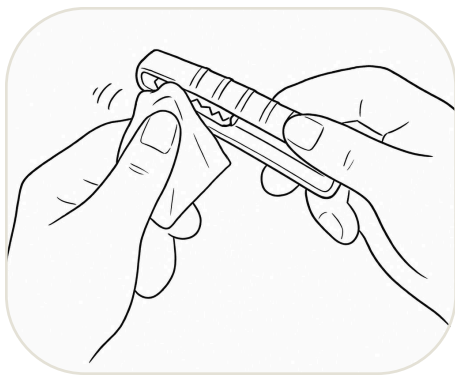


1 Clean Yewtwist

Wipe Yewtwist with a hospital-grade disinfectant wipe.




2 Select jaw size

Choose the jaw that best fits the connector. Place it on the connector furthest from the patient.

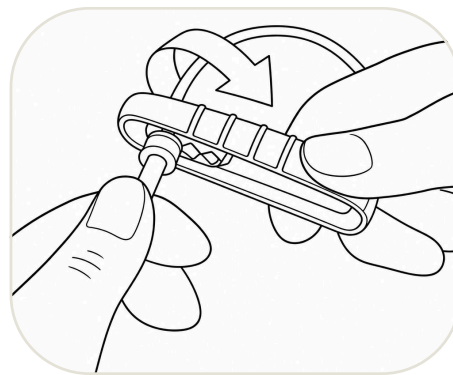

3 Pinch and turn

Pinch the arms together and turn counterclockwise.

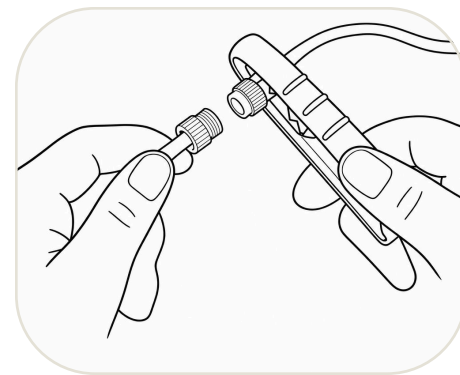

4 Detach

Separate the connection.

Cleaning instructions

Clean Yewtwist before and after each use with a hospital-grade disinfectant wipe. Most hospital-grade disinfectant wipes may be used.

- 1 Wipe all surfaces thoroughly, including the jaws and the grip ridges.
- 2 Follow the contact time in the wipe manufacturer's instructions for use.
- 3 Allow to air dry before use or storage.

Tested disinfectants

Diversey™ Oxivir® Tb RTU Disinfectant Cleaner

Clorox Healthcare® Bleach Germicidal Cleaner

Clorox Healthcare® Hydrogen Peroxide Cleaner

Super Sani-Cloth® Wipes

Metrex CaviWipes™

Inspection and replacement

Inspect Yewtwist before each use. Replace it if the jaws are cracked, torn, no longer grip, or if the body is cracked or bent.

Warnings, precautions and storage

- Nonsterile.
- Not made with natural rubber latex.
- Pinch point hazard.
- Do not apply excessive force when attempting to loosen a connection that is firmly seized or cannot be safely disconnected.
- Do not pull on the line or tube while Yewtwist is attached.
- Do not use Yewtwist to tighten connections. Connect by hand.
- Not a toy. Keep away from small children and pets.
- For hospital, clinic and home use.
- Store clean and dry at room temperature, in its bag.